



RiseAtlas450M™ with QuickTrolleySystem, high humidity SystemRoMedic™

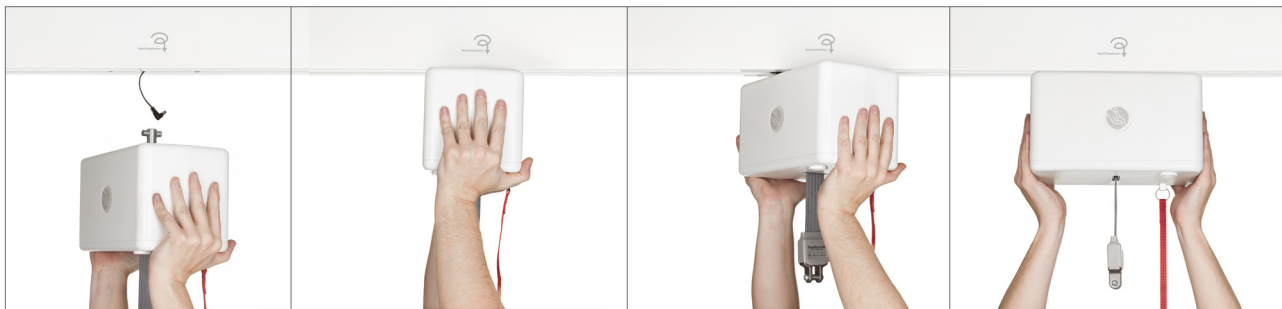

QuickTrolleySystem facilitates mounting and dismantling from the rail and provides the possibility to move the ceiling lift unit between different rail systems in different rooms for increased utilization and cost efficiency.

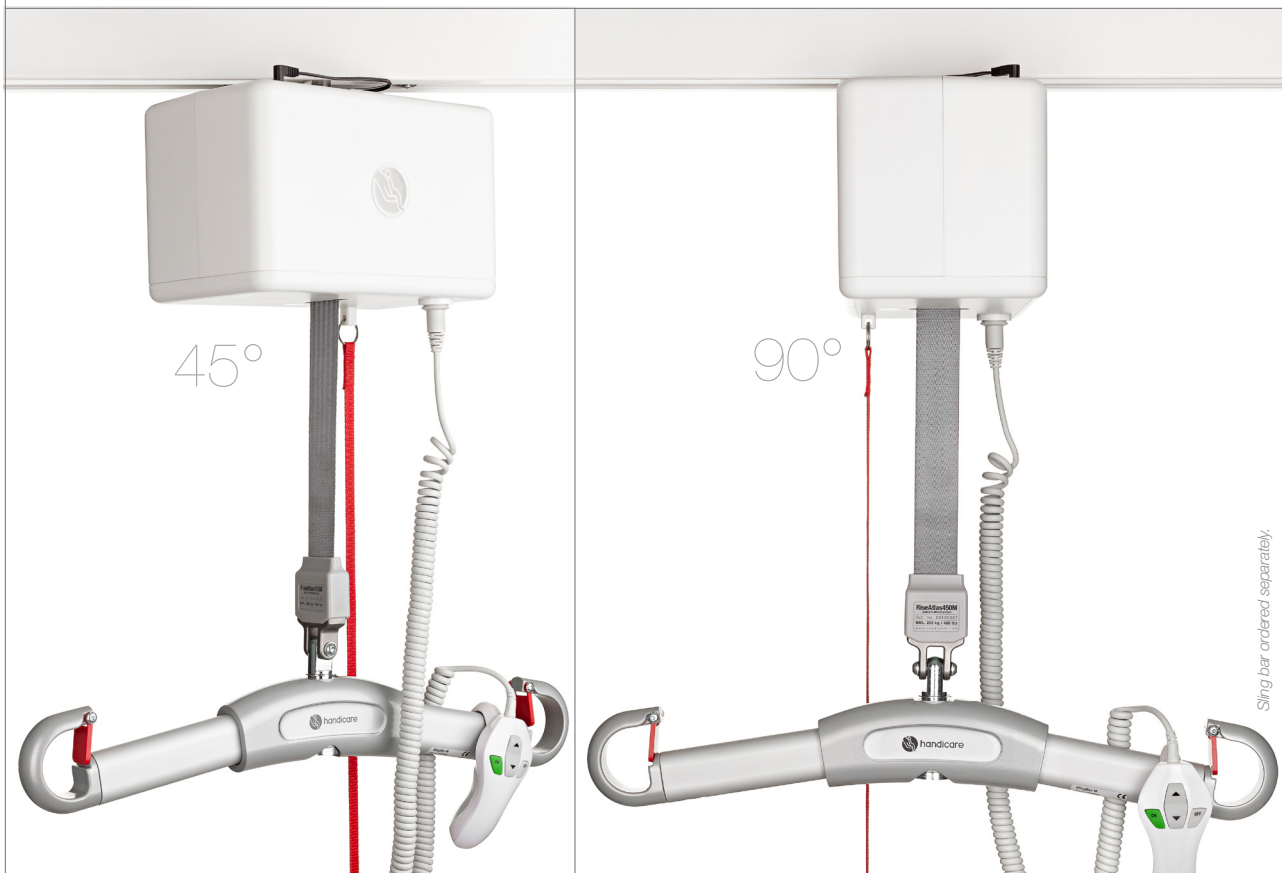

With QuickTrolleySystem, the ceiling lift unit and the lift strap can easily be preset to an angle of 45 or 90 degrees towards the rail for smooth and convenient transfers of users from one room to another.

Sling bar ordered separately.

Special version for wet environments – for users weighing up to 205 kg/450 lbs

- SWL 205 kg/450 lbs
- Manual transfer in rail
- Special version for use in environments where moisture and dampness are present
- Modern, elegant and hygienic design
- Unusually small measurements and extremely low weight
- For all rail systems, in all settings and in all types of rooms
- Large lifting area and possibility of both very low and very high lifting
- In all situations and for all user needs
- Handicare's unique SystemRoMedic™ QuickTrolleySystem provides more flexibility, versatility and convenience
- Very easy to use, providing safe working conditions and allowing the caregiver to remain close to the user throughout the entire lifting procedure
- Moves silently and gently, soft start and stop
- Ergonomic and easy-to-use hand control
- On/Off-button on the hand control prevents unauthorized use of lift functions
- Manual emergency stop
- Electrical and manual emergency lowering
- Overload protection
- Over speed protection
- Slack strap switch
- Night light
- Efficient strap control preventing accidental twisting of the lift strap
- High capacity, environmentally friendly NiMH battery
- Unique charging solution, with one trolley for two charging options (in-rail/end point charging)
- Discreet information panel
- Easy to service
- Diagnostics, a special software available as an accessory allows a PC to access further information about the ceiling lift unit
- Many accessories for safe, easy and individually adapted lifting of users in most transfer situations



RiseAtlas450M™ with QuickTrolleySystem, high humidity

For professional and highly functional care environments with a harmonic and human atmosphere.



Ergonomic and user-friendly hand control



PRODUCT SPECIFICATION

Item No. **Product**
 50100051 ___ RiseAtlas450M, with QuickTrolleySystem, high humidity
RiseAtlas450M, with QuickTrolleySystem, high humidity is supplied with hand control and charger. Sling bar ordered separately.

ACCESSORIES

Item No.	Product	Description
70200046 ___	QuickTrolley450	Quick release trolley, for 50100057, 50100051
70200051 ___	QuickTrolley450 70 mm	Quick release trolley, lowering lift unit by 70 mm / 2.76", for 50100057, 50100051
70200075 ___	QuickTrolleyUniversal 360	Quick release trolley, turns 360 degrees
70200048 ___	Software for diagnostics	Software, for all models and special versions of RiseAtlas
50400066 ___	Cover to charger	For all models and special versions of RiseAtlas
50400048 ___	Parking Placard	For convenient storing of sling bars, manuals and lifting slings
50400332 ___	Parking Placard BariBar	For convenient storing of BariBar and lifting slings
70200071 ___	SlingBar XS	250 mm / 9.8"
70200001 ___	SlingBar S	350 mm / 13.8"
70200002 ___	SlingBar M	450 mm / 17.7"
70200003 ___	SlingBar L	600 mm / 23.6"
70200042 ___	SlingBarSpreader M	Side bars for SlingBar XS/S/M/L
70200072 ___	BariBar	Four-point sling bar for ceiling lift units, SWL 454 kg/1000 lbs
70200017 ___	Sling bar RFL X4	Four-point sling crossbar, SWL 300 kg/660 lbs
70200053-___		
70200058 ___	Lift strap extension	Extension straps, 6 lengths: 370-970 mm / 14.6"-38.2"

TECHNICAL INFORMATION

Safe working load: _____ 205 kg / 450 lbs
 Material, chassis: _____ Aluminium
 Material, cover: _____ Non flammable ABS plastic
 Measurements, length x width x height: ___ 266 x 151 x 156 mm / 10.5 x 5.9 x 6.1"
 Weight: _____ 7.3 kg / 16.1 lbs
 Transfer in rail: _____ Manual
 Hand control: _____ Electrical, 4 buttons
 Lifting range: _____ 2400 mm / 94.5"
 Lifting speed, without load: _____ 50 mm / 2.0" per sec
 Lifting speed, with load: _____ 39 mm / 1.5" per sec
 Sound level, with load: _____ 55 dB(A)
 Efficient strap control: _____ Standard
 Emergency lowering: _____ Electrical and manual
 Emergency stop: _____ Manual
 Night light: _____ Standard
 Battery: _____ NiMH 24V 3.2 Ah/768Wh
 Hand control charging: _____ No
 In-rail charging: _____ Optional, standard
 End point charging: _____ Optional, standard
 Protection class, lift unit: _____ IPX4
 Protection class, hand control: _____ IPX4



RiseAtlas450M with QuickTrolleySystem, high humidity meets the requirements of the Medical Devices Directive for Class 1 products 93/42/EEC.

Please visit www.handicare.com for more information.

SystemRoMedic™ is a complete solution for easy transfers. The concept encompasses four product areas that meet most needs in different transfer situations.

Transfer, between two locations. **Positioning**, within the same location.

Support, for mobility. **Lifting**, both manual and mechanical.



handicare

Handicare AB
 Maskinvägen 17,
 SE 972 54 Luleå, SWEDEN
 Tel: 08-557 62 200
 Fax: 08-557 62 299
www.handicare.se
info@handicare.se