

MiniLift125

SystemRoMedic™

SWL
125 kg
275 lbs

Easily adjustable lift arm

Easy to manoeuvre

Grip-friendly handles and side supports provide multiple hand placement options for the user

Quick guide on the lift

Unique sit-to-stand pattern

High capacity battery with long service life

Compact format

Built-in charging

Soft, adjustable lower leg support


Lockable rear castors

Footplate with low step-in height and anti-slip surface

Electrical base widening

A unique solution for active sit-to-stand – with electrical base widening

MiniLift125 is a mobile sit-to-stand lift which is used, in combination with a specially adapted lifting vest, to assist users with reduced function when rising from a sitting to a standing position. MiniLift125 features a construction offering a totally unique and extremely comfortable sit-to-stand experience. The lift is also compact and convenient and very easy to manoeuvre; features which, in combination with the possibility to adjust the settings for an optimal fit for every user, makes it very useful in most locations and situations.

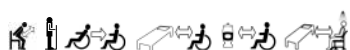

MiniLift125

Extremely comfortable, active sit-to-stand with natural movement pattern

MiniLift125 has been developed to, as gently and comfortably as possible, assist the user when rising from a sitting to a standing position. The handles and the side supports offer multiple hand placement options for the user. The stable anti-slip foot plate is very low and the lower leg supports, as well as the lifting arm, can be adjusted to fit every user optimally. When MiniLift125 is used in combination with the appropriate lifting accessories, the user is supported under the feet, for the front of the lower legs and behind the back, which provides for a safe and secure, yet active sit-to-stand procedure. The construction safely moves the user forwards and upwards in a natural movement pattern and, at the same time, the leg muscles and balance are exercised. MiniLift125 is manufactured in steel for stability and strength and is suitable for users weighing up to 125 kg/275 lbs.

Active sit-to-stand with natural movement pattern.



PRODUCT SPECIFICATION

MiniLift125 __ Mobile sit-to-stand lift with electrical base widening __ Art. No. 400641434

TECHNICAL INFORMATION

Safe working load (SWL)	125 kg / 275 lbs
User length, min-max	1400-2000 mm / 55.1-78.7"
Material	Powder-coated steel
Height	995 mm / 39.2"
Length	900 mm / 35.4"
Width of base, outer measurement	600-890 mm / 23.6-35.0"
Width of base, inner measurement	445-745 mm / 17.5-29.3"
Weight	39 kg / 86.0 lbs
Base widening	Electrical
Height of footplate, for user	50 mm / 2.0"
Castor size, diameter	80 mm / 3.1"
Lifting speed, without load	35 mm / 1.4" per sec
Emergency lowering	Manual and electrical



MiniLift125 meets the requirements of the Medical Devices Directive for Class 1 products 93/42/EEC.







Ergonomic and easy-to-use hand control with 4 buttons.




New control box with diagnostics, e.g. lift counter and overload indicator.

Emergency stop button readily available on the control box.

ACCESSORIES

	Art. No.	Product name	Description
 45500003-009	45500003	ThoraxSling	Polyester, X-Small
	45500004	ThoraxSling	Polyester, Small
	45500006	ThoraxSling	Polyester, Medium
	45500007	ThoraxSling	Polyester, Large
	45500008	ThoraxSling	Polyester, X-Large
 45550004-009	45500009	ThoraxSling	Polyester, XX-Large
	45550004	ThoraxSling	Wipeable, Small
	45550006	ThoraxSling	Wipeable, Medium
	45550007	ThoraxSling	Wipeable, Large
	45550008	ThoraxSling	Wipeable, X-Large
 45600004-008	45550009	ThoraxSling	Wipeable, XX-Large
	45590004	ThoraxSling	Disposable, non-woven, Small, 5 pcs
	45590006	ThoraxSling	Disposable, non-woven, Medium, 5 pcs
	45590007	ThoraxSling	Disposable, non-woven, Large, 5 pcs
	45590008	ThoraxSling	Disposable, non-woven, X-Large, 5 pcs
 45600004-008	45600004	ThoraxSling	With seat support, polyester, Small
	45600006	ThoraxSling	With seat support, polyester, Medium
	45600007	ThoraxSling	With seat support, polyester, Large
	45600008	ThoraxSling	With seat support, polyester, X-Large
	45690004	ThoraxSling	With seat support, disposable, non-woven, Small, 5 pcs
 45690004-008	45690006	ThoraxSling	With seat support, disposable, non-woven, Medium, 5 pcs
	45690007	ThoraxSling	With seat support, disposable, non-woven, Large, 5 pcs
	45690008	ThoraxSling	With seat support, disposable, non-woven, X-Large, 5 pcs
	70200033	CalfStrap	For MiniLift125 as from serial number 125001653
	770200089	Hand control HB33-6	6 buttons, with service, battery status and load indicators
	70200029	External charger Linak	External charger for all models and variants of MiniLift

**ThoraxSling**
SystemRoMedic™

ThoraxSling is a specially adapted and very comfortable sit-to-stand vest that is used in combination with Handicare's SystemRoMedic™ MiniLift sit-to-stand lifts. It is easy for the caregiver to apply, while giving the user an added sense of security during rising and standing training. ThoraxSling also facilitates toileting. ThoraxSling is available in two models; one which provides support behind the back (ThoraxSling) and one that also gives support under the seat (ThoraxSling, with seat support), and in several materials and sizes.

ThoraxSling

SystemRoMedic™

ThoraxSling is a specially adapted and very comfortable sit-to-stand vest that is used in combination with Handicare's SystemRoMedic™ MiniLift sit-to-stand lifts. It is easy for the caregiver to apply, while giving the user an added sense of security during rising and standing training. ThoraxSling also facilitates toileting. ThoraxSling is available in two models; one which provides support behind the back (ThoraxSling) and one that also gives support under the seat (ThoraxSling, with seat support), and in several materials and sizes.

Other models and variants of MiniLift mobile sit-to-stand lifts

MiniLift is available in several models and variants with different lifting capacity, with electrical or manual base-width adjustment and with standard height or low legs for use with lower beds.

Art. No.	Product name	Description
401100334	MiniLift200	Mobile sit-to-stand lift with electrical base widening, SWL 200 kg/440 lbs
60300010	MiniLift160EM	Mobile sit-to-stand lift with manual base widening, SWL 160 kg/350 lbs
60300011	MiniLift160EM	Mobile sit-to-stand lift with manual base widening, low legs, SWL 160 kg/350 lbs
60300012	MiniLift160EE	Mobile sit-to-stand lift with electrical base widening, SWL 160 kg/350 lbs
60300013	MiniLift160EE	Mobile sit-to-stand lift with electrical base widening, low legs, SWL 160 kg/350 lbs

Please visit www.handicare.com for more information.



handicare

Maskinvägen 17,
972 54 Luleå, Sweden
Tel: +46 (0)8-557 62 200
Fax: +46 (0)8-557 62 299
www.handicare.com
info@handicare.se

SystemRoMedic™ is a complete solution for easy transfers. The concept encompasses four product areas that meet most needs in different transfer situations.

Transfer, between two locations. **Positioning**, within the same location.

Support, for mobility. **Lifting**, both manual and mechanical.